PERSONENWAAGE

Art.-Nr.: 74999 / ZFC5020

DE / EN / NL

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Size unit: cm, feet, inches, inches

Memory: 10 users

User modes: adult, child (not suitable for children under 6 years old)

Recommended operating ambient temperature: 0- 40 deg.

Celsius /32 - 104 degrees Fahrenheit

Relative Humidity - 85%

INSTRUCTIONS

This scale is for home use only, not for commercial use.

Inserting the Battery

- Open the battery compartment under the scale and insert the batteries.
- Insert the batteries according to the battery polarity. Close the battery compartment.


Selecting the unit of measurement

- Please select the unit of measurement for the weight - you can choose between kg and lb.

First steps

- Insert the battery into the battery compartment according to the polarity.
- Turn the scale on by pressing the key (see diagram on the next page), and use the up/down keys to select a unit of measure.
- Place the scale on a stable, level surface to ensure accuracy.
- Step on the base of the scale to turn it on.
- Stand still while the scale calculates your weight.
- Your weight will be displayed and flash twice before it is permanently displayed.
- When you step on the scale, the scale will automatically turn off after 8 seconds.
- If no weight is recorded, 0.0 is displayed and The scale will automatically turn off after 8 seconds.

OPERATION

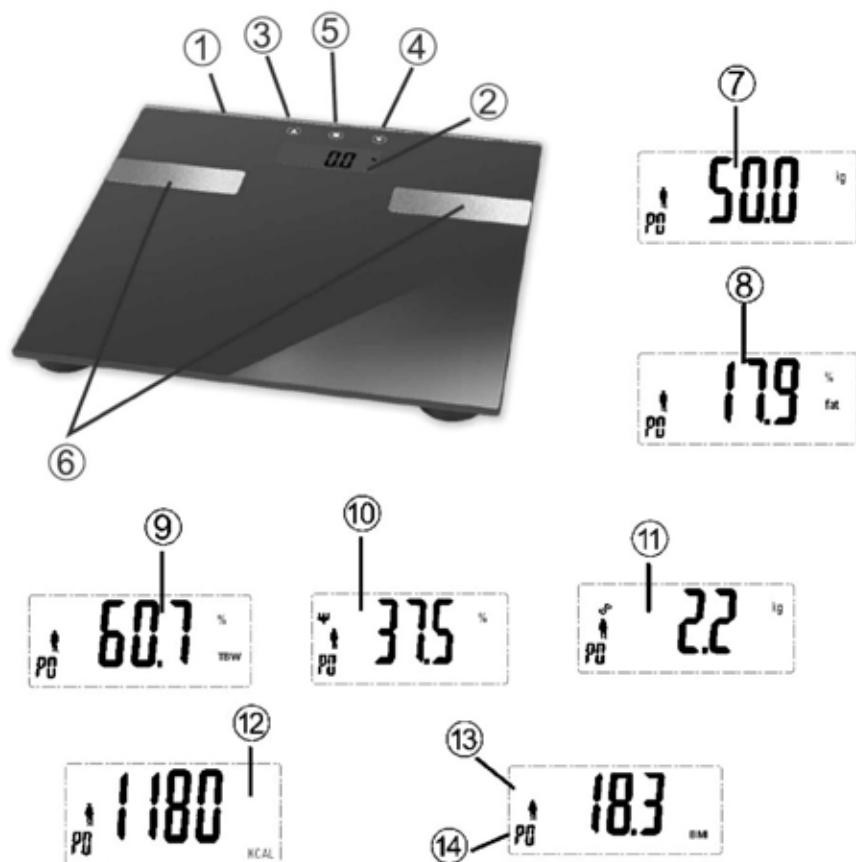1. battery compartment (on the bottom side)
2. display
3. UP button ▲ (on top, under the glass platform)
4. DOWN button ▼ (on top, under the glass platform)
5. SET button (on the glass stand - to store your personal data)
6. electrodes (you have to step on them barefoot to check your body fat)
7. body weight value
8. body fat value
9. body water value
10. muscle mass value
11. bone mass value
12. calorie requirement value
13. BMI value
14. default setting P0 (e.g. user profile 0, female)

WARNING: DANGER OF SLIPPING!

Danger of slipping if feet are wet!

ATTENTION!

- Gently step on your scale and the weight will be displayed automatically. If the scale is not used, it will turn off automatically.
- Remove the battery if the scale will not be used for an extended period of time.
- Clean the scale with a soft cloth moistened with water or a mild detergent. Never clean the scale with gasoline, hot water or harsh chemicals.
- Do not expose the scale to extreme temperatures, humidity, direct sunlight, moisture or dust.
- The batteries supplied with the scale are for demonstration purposes only. We recommend that new batteries be used.
- When the battery level is low, the „Lo“ symbol is displayed.
- Do not subject the scale to shock or impact, drop it to the floor, or attempt to open it and tamper with the mechanisms or circuitry.
- Our scale is designed for personal home use only. It is not suitable for professional use in hospitals, medical facilities or retail.

SPECIFICATIONS

Auto-off: If the scale is not used for about eight seconds, it will turn off automatically.

- Power supply: 1×3V CR2032 battery.
- Unit: kg, lb, st
- Weighing range: 5 kg - 180 kg
- Low battery indicator: Lo
- Overload indicator: Err
- Graduation: 100g
-

WARNING:

If the battery belonging to the item is swallowed, this may result in serious internal injuries and even death.

- Keep both new and used button cell batteries out of reach of children.
- If the battery compartment can no longer be closed, stop using the item and store it out of the reach of children.
- If you suspect that a button cell battery has been swallowed or otherwise entered your body, seek medical attention immediately.

When moving the scale, perform the following initialization:


Step 1 - Press the center of the platform by touching it with your hand or foot.

Step 2 - "0.0" will be displayed.

Step 3 - The scale will turn off and is now ready for use.

Step 4 - Carefully step on the scale and stand on it without moving. Wait while the scale calculates your weight.

The correct weight will then be displayed as follows: 80.3 kg.

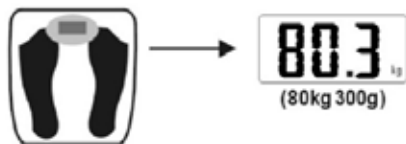WEIGHING AND MEASURING

To ensure accurate results, the user should enter the required personal data.

Parameters: Age, gender and height.

This scale can store the personal parameters of 10 people.

The Bioelectric Impedance Analysis (BIA) technology used in this scale is more technologically advanced than older Body Mass Index (BMI) methods that calculate fat value from height and weight only.

The BIA technology uses 5 parameters to provide a more accurate resolution for the analysis (weight and the biological resistance and the 3 personal, manually entered parameters).

BODY FAT MEASUREMENT:

Mode selection and parameter setting:

Press ■ to turn the scale on or off.

Press ▲ or ▼ to switch between body parameters 0-9. Press to confirm ■

Press ▲ or ▼ to select a man or a woman. Press to confirm ■

Press ▲ or ▼ to specify the size. Press to confirm ■

Press ▲ or ▼ to specify the age. Press to confirm ■


The parameter settings flash and are recorded. Step on the scale barefoot when "0.0" appears.

Place both feet firmly on the electrodes.

Weight is measured, recorded and displayed. When complete, weight and body fat percentage, water percentage, bone weight, muscle percentage, KCAL and BMI are displayed three times alternately before the scale automatically turns off.

The body fat measurement on the opposite page (body fat and water percentage table) provides information about your body condition.

The screen sequence is as follows:

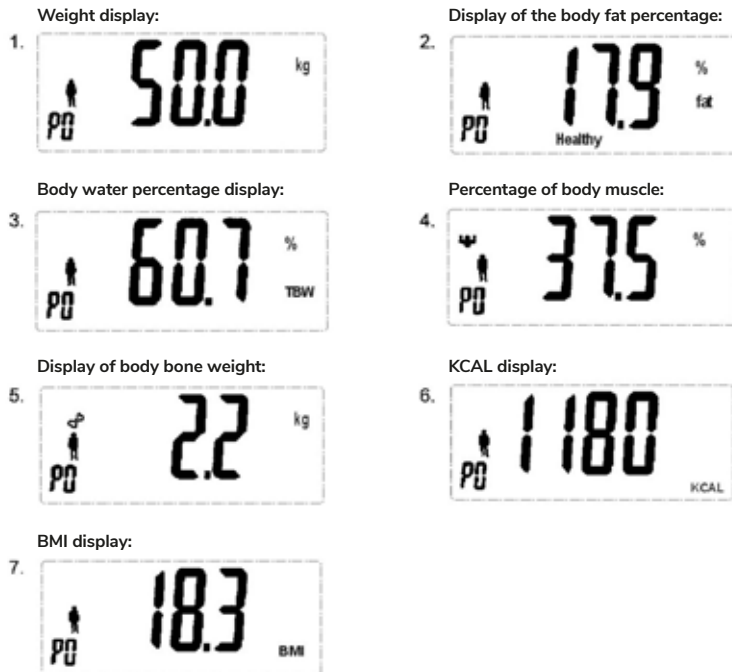ERROR DISPLAY

- If "LO" appears on the LCD, you must install a new battery.
- If the LCD displays "ERR", the maximum weight has been exceeded.
- If the LCD displays Err 2, there is an error in the body fat percentage.

ATTENTION:

- Weigh yourself without clothes and shoes, before meals and always at the same time of day.
- For better accuracy, do not weigh yourself within 2 hours of waking up.
- Keep the scale clean and dry to protect electronic components.
- Clean with a damp cloth after each use, do not use solvents or immerse in water.
- Do not store the scale in an upright position when not in use or the battery may leak.
- Avoid shocks or vibrations to the device, e.g. by dropping the scale on the floor. Do not place any objects on the scale.
- A user who has metal parts implanted in his body may get different results when measuring body fat.
- Wearers of pacemakers should not use the device.

DISPOSAL

The device (including removable parts and accessories) must not be disposed of with the residual waste at the end of its useful life, but must be disposed of in accordance with the European Directive. Since the scale must be disposed of separately from household waste, either take it to an authorized disposal point for electrical and electronic equipment or return it to the dealer when you purchase a new instrument for the same purpose. Any violation will be strictly punished. The batteries used in the device must be disposed of in special containers when they are used up.

 **NOTE:** Specifications and designs are based on the latest information. available at the time of available at the time of printing and are subject to change without notice.

REFERENCE TABLE

Table body fat ratio (%)

AGE	Female (%)				Male (%)			
	Under-weight	Healthy	Over-weight	Obese	Under-weight	Healthy	Over-weight	Obese
6-39	<21	21≈33	33≈39	>39	<8	8≈19	19≈25	>25
40-59	<23	23≈35	35≈40	>40	<11	11≈22	22≈28	>28
60-100	<24	24≈36	36≈42	>42	<13	13≈25	25≈30	>30

Due to a natural body variation between people, this table may only be used for reference purposes.

Water Ratio Table & Muscle Ratio Table (%)

AGE	Female (%)				Male (%)			
	Low	Healthy	High	Body Muscles	Low	Healthy	High	Body Muscles
6-15	<21	57,0-67,0	>67,0	>34	<58,0	58,0-72,0	>72,0	>40
16-30	<23	47,0-57,0	>57,0		<53,0	53,0-67,0	>67,0	
31-60	<24	42,0-52,0	>52,0		<47,0	47,0-61,0	>61,0	
61-100	<24	37,0-47,0	>47,0		<42,0	42,0-56,0	>56,0	

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Hatex AS GmbH & Co. KG