MICHELINO®
— since 1951 —

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KÖRPER ANALYSE WAAGE

Art.-Nr.: 75050 / 75051

DE / EN / FR / NL

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E-Mail: info@hatex24.de

IMPORTANT INFORMATION BEFORE USE

Read the operating instructions and safety instructions carefully before using the device!

- Failure to observe the operating instructions may result in serious injury or damage to the device.
- Keep the operating instructions for further use.
- If you pass the device on to third parties, be sure to include these operating instructions.
- The device is designed for private household use only and not for commercial use.
- Only use the device indoors and not outdoors.

WARNINGS

The following warnings are used in this manual:

 **DANGER! High risk:**

Disregarding this warning can cause harm to life and limb.

 **WARNING! Average risk:**

Disregarding this warning can cause injury or serious property damage.

 **CAUTION: Low risk:**

Disregarding this warning may cause minor injury or property damage.

NOTE: Facts and special features that should be observed when handling the device.

SAFETY INSTRUCTIONS

- If the safety instructions are ignored, the manufacturer cannot be held liable for any damage.
- Keep these instructions in a safe place for future reference.
- Always place the scale on a stable and level surface to ensure correct measurement. If the surface yields (carpet, etc.), the four force transducers in the feet will be unevenly loaded and the measured values will inevitably be falsified.
- Do not use abrasive cleaning agents for cleaning.
- Remove the batteries if the scale will not be used for a long period of time. Store the scale in a warm and dry place.
- Do not attempt to dismantle or repair the scale (loss of warranty!).
- The scale contains sensitive electronic components. Keep the scale away from intense heat and environments with high humidity.

- Always handle the scale with care. (Never throw or jump on it)
- The scale is designed for home use only and for measuring people.

⚠ WARNING!

- Do not step on the scale with wet feet –  **danger of slipping!**

⚠ WARNING!

- This product is not suitable for people wearing a pacemaker or another internal medical device.

⚠ CAUTION!

- This product should not be used to diagnose or treat any medical condition.

BATTERY

1. The battery compartment is located on the underside of the unit.
2. Open the lid and insert the batteries into the battery compartment.
3. Close the lid of the battery compartment again.

⚠ DANGER!

Swallowing the battery that comes with the article can cause serious internal injury and even death.

Keep both new and used batteries out of reach of children.

If the battery compartment can no longer be closed, do not use the article and keep it out of the reach of children.

If you suspect that a battery has been swallowed or otherwise entered your body, consult a doctor immediately.


NOTE:

When a battery symbol appears in the upper left-hand corner of the display, the inserted battery must be replaced with a full one.

SCOPE OF DELIVERY

1. Body Analysis Scales
2. Instruction manual

DEVICE DESCRIPTION

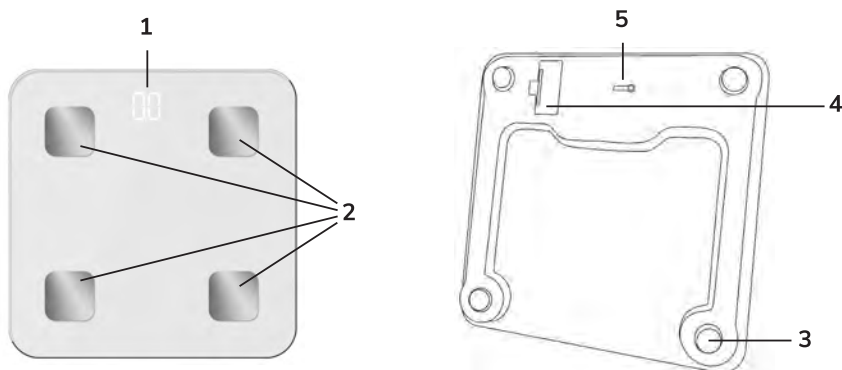1. LED-Display
2. Biometric Sensors
3. Non-slip feet
4. Battery compartment
5. Switch for setting the "kg" / "lb" unit

TECHNICAL DATA

Model:	28406 / 28407 / 28408
Battery:	2*AAA*1,5V
Operating temperature:	5–35 °C
Maximum load capacity / weighing range:	180 kg / 396 lb
Graduation:	1 kg / 0.2 lb
Measuring unit:	"kg" / "lb"

BEFORE FIRST USE


⚠ CAUTION!

Packaging material is not a child's toy. Children must not play with the plastic bags. There is a risk of suffocation.

- Remove all parts from the packaging and remove all packaging materials. Also remove any transport locks.
- Check the scope of delivery for completeness and damage and clean the unit.

COMMISSIONING

1. Visit designhq.com/app-download to download the iHome Health & Fitness app. You can also scan the QR Code:

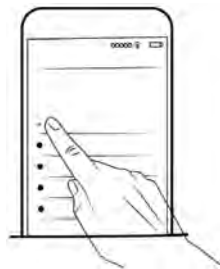2. Pair the iHome Smart Scale + with the iHome Health & Fitness App. Once your account is created, you can head to the devices tab. The devices tab will be the 3rd tab on the bottom of the app. From here, tap "Add Device" and look for the iHome Smart Scale +.
3. When stepping on the Smart Scale +, ensure the app is open to sync your weight and other body metrics. To ensure an accurate reading, ensure you step on the scale and cover the sensors with your bare feet. **Please note:** The scale must be placed on a flat, hard surface. If it is not, the readings may fluctuate.
4. The iHome Smart Scale + first measures your weight. If you continue to step on the scale with your bare feet, the iHome Smart Scale will begin to show dots on the LED. This indicates that it is currently measuring your body metrics.

⚠ CAUTION!

- Do not place heavy objects on the scale for a long time when you are not using the scale.

MEASUREMENTS

Body fat weight

Body muscle percentage

Body protein percentage

Weight

Health age of body

Obesity level

Weight of fat around organs

Body fat percent

Body health

Body water percentage

Water weight

Weight of bone

Body weight without fat

Jiggly fat weight

More...

INSERT OR REPLACE BATTERY

- When the message “Lo” appears in the display, the inserted battery must be replaced with a full one.
- The battery compartment is located on the underside of the machine.
- Open the lid and insert the batteries into the battery compartment.
- Close the lid of the battery compartment again.

CLEANING

- Use only a cloth dampened with water to clean your scale.
- The use of any detergents or cleaners may cause irreversible damage to the coating that protects the electrodes (2) used to measure your body fat.

STORAGE

- Store the device in a clean, dust-free and dry place

REQUIREMENTS

FC This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

DISPOSAL

Notes on environmental protection

The symbol opposite indicates that this product is subject to Directive 2012/19/EU. This directive states that you must not discard this product at the end of its working life with your normal household waste at the end of its service life, but must be collection points, recycling centres or waste management disposal companies.

Batteries and rechargeable batteries must not be disposed of with household waste!

As a consumer, you are legally obliged to dispose of all batteries and accumulators, regardless of whether or not they contain harmful substances*, or not, to a collection point in you municipality/borough or to a retailer, so that they can be disposed of in an environmentally friendly manner.

*marked with: Cd = cadmium, Hg = mercury, Pb = lead

For addresses and opening hours of collection points your local administration will be pleased to inform you.

WARRANTY

You receive a 3-year guarantee on this device from the date of purchase.

WARRANTY CONDITIONS

- The guarantee period begins on the date of purchase. Keep the original sales receipt as proof of purchase.
- Should a material or production defect occur within the three years from the date of purchase, we will, at our option, replace the product or repair it free of charge, provided that both the device and the receipt, including a brief written description of the defect, are submitted to us within the aforementioned three-year period.
- If the defect is covered by our guarantee, you will receive a new product or we will repair the defective device.
- The replacement or repair of the product does not start a new guarantee period.
- Replacing or repairing parts of the device does not extend the guarantee period. The warranty period is not extended by the guarantee.
- Damage or defects already present at the time of purchase must be reported immediately after unpacking the goods. Repairs after the expiration of the guarantee period are subject to a charge.
- Covered by the guarantee period is the product and the receipt, including a brief written description of the defect. Also, easily breakable or easily damaged components, such as batteries or parts made of glass, are not covered by the guarantee.
- If the defect is covered by our guarantee, you will receive a new product or we will repair the defective product. Take and follow all information and instructions for proper use and maintenance of the product from the operating instructions.
- The guarantee period is not extended by the guarantee.
- Damage or defects already present at the time of purchase must be reported immediately after unpacking the goods.
- Repairs after the expiration of the guarantee period are subject to a charge.
- Covered by the warranty are material or manufacturing defects.
- Not covered by the guarantee are parts of the product that are subject to normal wear and tear and are considered wear parts. Also, easily breakable or easily damaged components, such as batteries or parts made of glass, are not covered by the guarantee.
- The guarantee is void due to external damage to the product, improper use and improper maintenance. Take and follow all information and instructions for proper use and maintenance of the product from the operating instructions.
- The guarantee also does not cover defects or damage caused by failure to follow the operating instructions.

- Uses and actions advised against or warned against in the operating instructions must be avoided at all costs.

PROCESSING IN THE CASE OF A WARRANTY CLAIM

To ensure that your claim is processed quickly, please follow the instructions below:

- Please refer to the type plate of the device, the sales box or the operating instructions for the article number.
- If functional errors or other defects occur, please contact the following service department by telephone or e-mail. There, the further procedure will be described to you in detail.

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Hatex AS GmbH & Co. KG